








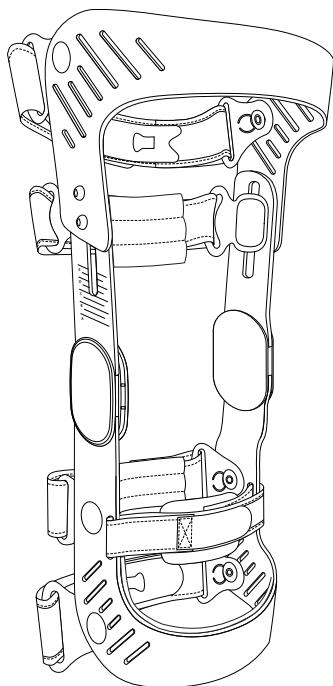
Instructions for Use

REBOUND[®] DUAL

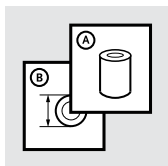		3
EN	Instructions for Use	5
DE	Gebrauchsanweisung	9
FR	Notice d'utilisation	14
ES	Instrucciones para el uso	18
IT	Istruzioni per l'uso	23
NO	Bruksanvisning	28
DA	Brugsanvisning	32
SV	Bruksanvisning	36
EL	Οδηγίες Χρήσης	40
FI	Käyttöohjeet	45
NL	Gebruiksaanwijzing	49
PT	Instruções de Utilização	53
PL	Instrukcja użytkowania	58
CS	Návod k použití	63
TR	Kullanım Talimatları	67
RU	Инструкция по использованию	71
JA	取扱説明書	76
ZH	中文说明书	80
KO	사용 설명서	84

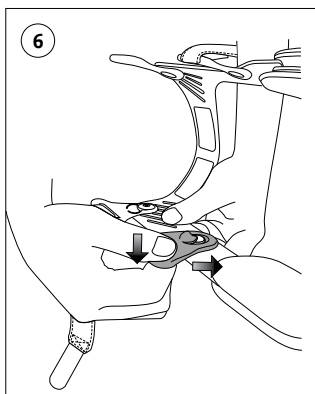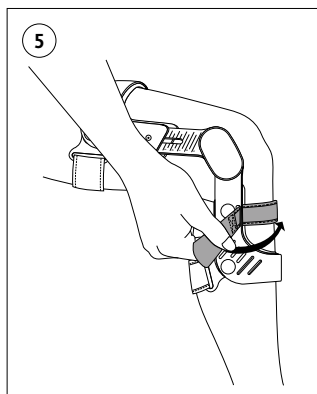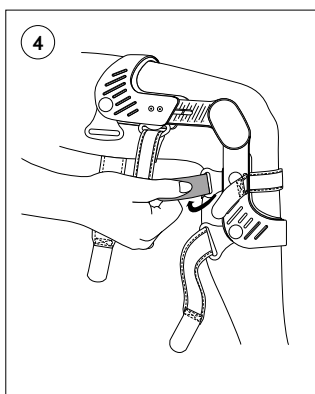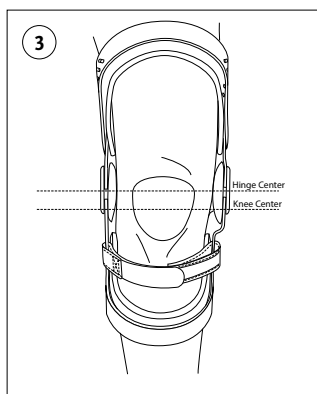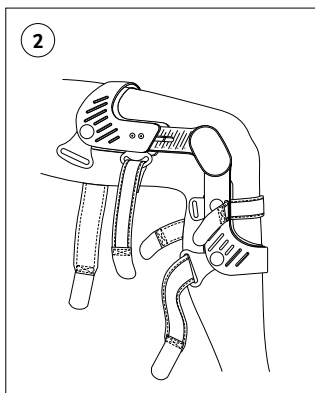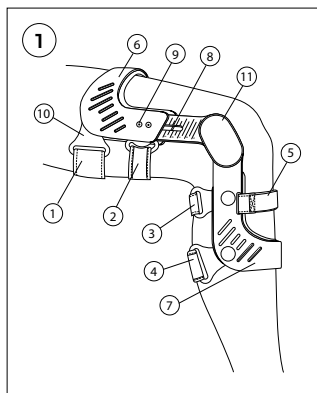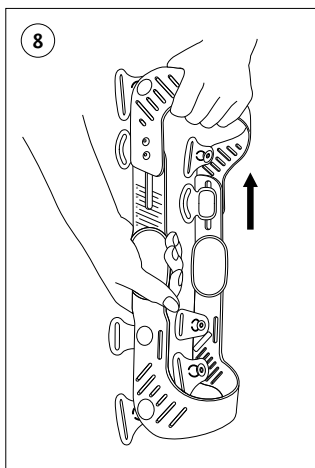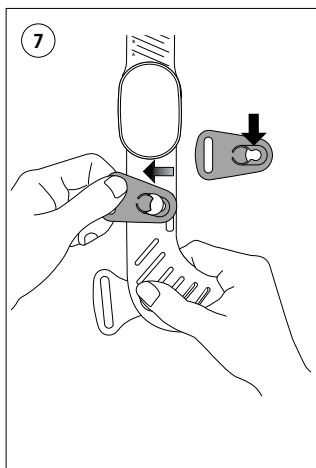ENGLISH

LIST OF PARTS: (Figure 1)

1. Top Strap
2. Upper Cruciate Strap
3. Lower Cruciate Strap
4. Bottom Strap
5. Anterior Tibial Strap
6. Upper Frame
7. Lower Frame
8. Acculign™ Hinge Arm
9. Acculign™ Hinge Arm Screw
10. D-ring
11. Condyle pad

INDICATION FOR USE

- Anterior Cruciate ligament instability, posterior cruciate ligament instability, medial collateral ligament injury, lateral collateral ligament injury.
- Conditions requiring unicompartamental load control.

WARNING

This knee brace is a supportive brace only and is not intended or guaranteed to prevent knee injury.

Purchaser agrees to the following:

1. Össur accepts no responsibility for any injury which might occur while using Össur knee bracing products.
2. With use of this type of device, it is possible to receive various types of skin irritation, including rashes and pressure points. If you experience any unusual reactions while using this product, consult your health care provider immediately.
3. If you have any questions as to how this knee brace applies to your own specific knee condition, consult your health care provider.

If you do not agree to the above conditions, contact your brace provider before using this product.

BRACE APPLICATION

1. Unfasten all straps of the brace.
2. With the knee bent in approximately 90 degrees, place the brace over your knee such that the patella is centered between the condyle pads. **(Figure 2)**
Note: This places the hinge slightly above the joint line so the brace can properly seat itself when you stand. **(Figure 3)**
Note: All straps can be trimmed or adjusted to length.
3. Fasten the lower cruciate strap just below the knee. Fasten the strap snugly above the calf muscle as this will aid in the brace suspension. **(Figure 4)**
4. Snugly fasten the bottom strap by inserting it through the D-ring.
5. Snugly fasten the upper cruciate strap above the knee.
6. Snugly fasten the top strap. **Note:** Over-tightening the two top straps can cause the brace to migrate.
7. Snugly fasten the anterior tibial strap **(Figure 5)**. If the brace is too anterior, loosen the anterior tibial strap and tighten the lower cruciate and bottom straps. If the frame is too posterior, loosen the lower

cruciate and bottom straps AND tighten the anterior tibial strap. Proper anterior/posterior placement will help prevent migration and ensure proper brace suspension.

Warning: If brace causes discomfort or irritation, discontinue use immediately and contact your health care provider.

BRACE MAINTENANCE

The Rebound DUAL brace is designed for minimum maintenance and care. Follow these easy maintenance steps to keep your brace in good condition.

CLEANING

- Liners, straps and pads may be removed for cleaning
- Do NOT machine wash
- Do NOT place in dryer
- Hand-wash using mild non-detergent soap
- Rinse thoroughly
- Air-dry

If migration occurs, clean lower liners and pads with an alcohol wipe or damp cloth. If brace is used in salt water or chlorinated water, rinse well and air dry.

HINGE MAINTENANCE

- Do NOT disassemble the hinge
- Check regularly for foreign materials in hinge (i.e. dirt or grass)
- Carefully clean dirt from hinge using water or compressed air
- Lubricate with silicone spray or dry graphite if needed

ACCESSORIES

The following accessories are available for your Rebound DUAL brace at an additional charge.

To order call your Össur Customer Service Department.

- SmartDosing kit
- Flexion stop kit
- Liner kit
- Strap kit
- D-ring kit
- Extension Stop Kit

SETTING UP THE D-RING

The field-serviceable d-rings can be changed from internal to external d-rings or replaced if a d-ring is broken. The following instructions show reconfiguration of internal to external d-rings.

1. Remove strap from d-ring.
2. Bend d-ring and push inward on internally mounted d-ring to disengage keyhole lock to remove from rivet. **(Figure 6)**
3. Replace d-ring to external configuration by pushing down onto rivet through the keyhole lock and sliding to click the d-ring into place. Keyhole lock should be placed with indentation side face up. **(Figure 7)**
4. Replace strap.
5. Repeat on other d-rings as needed.

REFURBISHING

From time to time, your Rebound DUAL should be sent to Össur for refurbishing. The frequency of this will depend on how often the brace is worn and for what activities. Your brace supplier should help you to determine when your brace needs to be sent in for service, but once a year is a good rule of thumb.

As part of the refurbishment process, Össur will check your brace for signs of wear and rebuild it with all new components, including rivets, buckles, straps and padding. There is a usual turnaround time of three to five working days. If you have any questions, contact your local Customer Service Department.

FIT ADJUSTMENTS

Contouring the Brace:

The Rebound DUAL brace can be adjusted for individual anatomies. The aircraft grade aluminum frame can be cold formed to accommodate different thigh and tibia shapes using bending irons as needed.

Adjusting the Varus / Valgus angle of the Brace:

The Rebound DUAL brace can be adjusted for individual varus or valgus angles and to create an offloading force to either the medial or lateral compartment. To adjust the angle, loosen the screws on the desired side of the brace and slide the upper frame to the desired position with a steady force. Adjust the medial upper frame to offload the medial compartment or adjust the lateral upper frame to offload the lateral compartment. The slide bar has a lettered scale and may be moved to increase or decrease the amount of varus or valgus force on the leg. **(Figure 8)** The lettered scale does not relate to degrees of offloading and should be used as a guide. Retighten the screws until snug (recommended torque value of 3 in-lb). All adjustments should be determined by a healthcare provider, in accordance with the patient's comfort level. Each time the brace is adjusted, the healthcare provider should note the letter for a record of progression.

Adjusting the Height of the Brace:

The Rebound DUAL brace can be adjusted for individual height requirements. The brace height can be adjusted by loosening the screws on both sides of the upper frame and pulling the frame up with a steady force. Ensure the brace is properly aligned by checking that both hinge arms are located on the same position on the medial and lateral hinge arm. Retighten the screws until snug.

Widening or narrowing fit at knee joint:

The Rebound DUAL brace can be adjusted for individual knee width. To widen or narrow the width at the knee, the medium condyle pads can be removed and replaced with the thick or thin version shipped with the Rebound DUAL. Any combination of pads can be used to achieve a more comfortable fit.

Changing Extension Stops

Extension stops are available for varying degrees of extension control.

1. Remove the pad from the condylar area.
2. Use a short Philips head screwdriver and remove the retaining screw.
3. Slide the 0° stop out of the hinge.
4. Insert desired extension stop.
5. Secure with retaining screw (DO NOT OVER TIGHTEN)

6. Repeat on the other hinge.
7. To add flexion stop: Follow step 1, insert desired flexion stop, and then follow steps 5–6.

Attaching PCL Strap

1. Note height position on medial and lateral hinge arms.
2. Remove hook-and-loop fastener to access connection area.
3. Remove Acculign hinge screws.
4. Remove Acculign hinge d-Ring.
5. Position PCL Acculign hinge d-ring.
6. Position PCL Acculign hinge inserts with raised portion into recessed pocket.
7. Return frame to initial positions on hinge arms.
8. Tighten Acculign screws.
9. Reapply hook-and-loop fastener.
10. Attach PCL strap.

LIMITED WARRANTY

This warranty covers all defects in materials and workmanship of your new Rebound DUAL Knee Brace. This warranty covers the frame and hinges for 2 years, and custom warranty is lifetime. The warranty covers all other materials, including straps, sleeves, liners and pads for six months from the date of original purchase.

Össur will repair or replace (at its option) any defective part covered under warranty at no charge. Repairs necessitated by accident, negligence, misuse or abuse by the purchaser are not covered. Also, consequential and incidental damages are not recoverable under this warranty. These limitations may not apply as some states do not allow the exclusion or limitation of incidental or consequential damages. Only the original purchaser is entitled to claim service under this warranty.

To obtain service under this warranty call our Customer Service Department.

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